



Specimen Inquiry for Beaker AP Users

1. Click on **Specimen Inquiry by Patient** from your Dashboard:

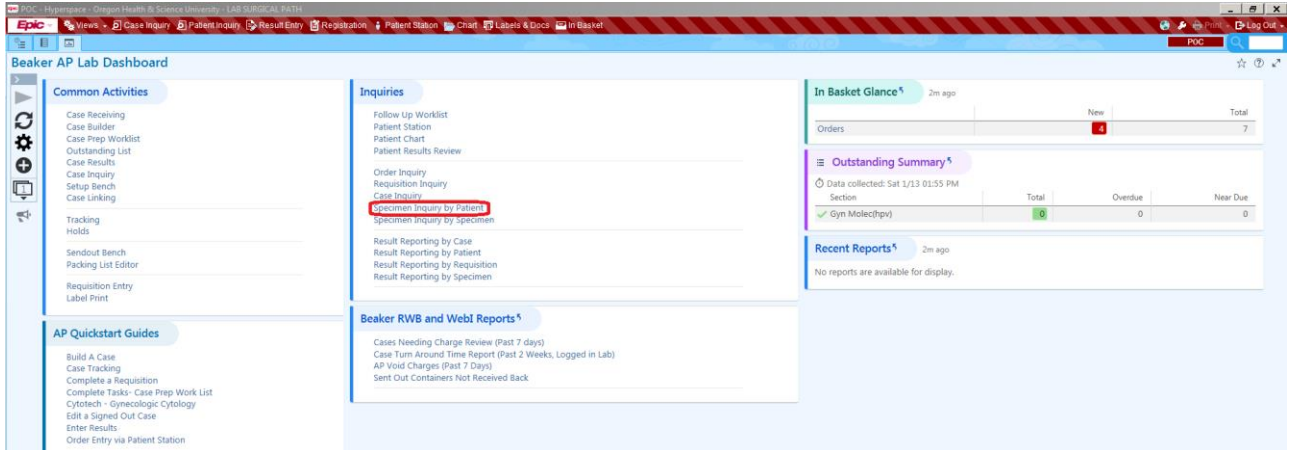


2. Search for the patient using name or MRN:

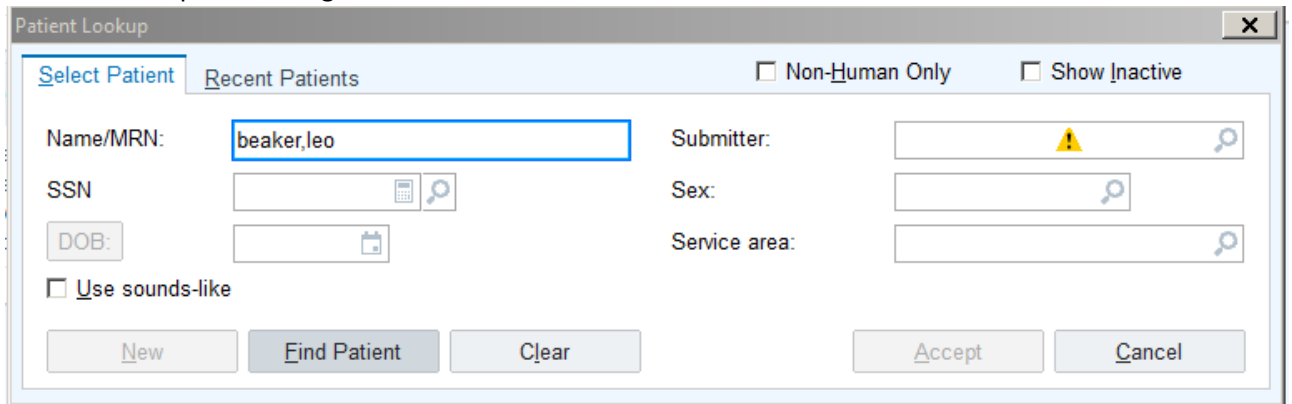

3. You will see a list of all the specimens for the patient with newest at the top. The tests done on each specimen will show to the right. The symbol under the "S" column will indicate the status of the test.

Specimens (Count: 23)	Specimen	Test	F	P	I	S	F/U	Submitter/Department	Requisition
1/8/2018 1420	*18KD-008P0001	GENETRAILS SOLID TUMOR PANEL		█	!	█		PROVIDENCE PORTLAND	*RQ6038
10/24/2017 1300	*17LP-297L0002	Plasma Lipids - Lipid Lab		█		×		ADVENTIST MEDICAL CENTER LABORATORY	*RQ4585
10/10/2017 1309	*17RF-283R0001	Miscellaneous Referral Test		█		×		LAB PHLEB PPV	
10/3/2017 1651	*17IM-276T0002	Lupus Inhibitor Evaluation Heparin, Either STD/Lmw APTT 1:1 Mix		█	!	█		PEACE HEALTH LABORATORY	*RQ4227
9/29/2017 1329	*17CL-272G0002	Fibrinogen Level		█		×		LAB PHLEB PPV	
9/11/2017 0910	*17ML-254M0001	Culture, Urine Susceptibility-Gram Negative MIC		█	!	█		MCMC Generic Submitter	*RQ3694
7/10/2017 1206	*17ML-191M0001	Culture, Urine Susceptibility-Gram Negative MIC Susceptibility-Gram Negative MIC		█	!	█		MCMC Generic Submitter	*RQ2309

- To see details about any specimen, click on the blue specimen number to open it up:

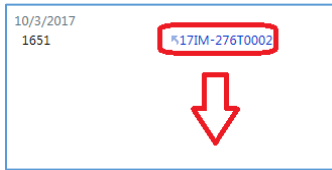

171M-276T0002 Instrument ID: 16560

Beaker, Leo (MRN 03427183) **Blood, Blood**

M. 14 yrs. 9/1/2003 Collected 10/3/2017 1651
 %RQ4227 submitted by PEACE HEALTH LABORATORY Container: 3 LT BLUE-2.7mL

Specimen Tracking

171M-276T0002.1 LT BLUE-2.7mL

Date	Detail	User	Location
10/03/2017			
1700	Result Filed To Chart	Service Account, Beaker Results Out	LAB SIC
1700	Charge Triggered	Mary Ulmschneider	LAB SIC
1700	Result Final Verified	Mary Ulmschneider	LAB SIC
1700	Test Canceled	Mary Ulmschneider	LAB SIC
1700	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1700	Result Entered	Mary Ulmschneider	LAB SIC
1700	Result Entered	Mary Ulmschneider	LAB SIC
1659	Charge Triggered	Mary Ulmschneider	LAB CORE
1659	Result Filed To Chart	Service Account, Beaker Results Out	LAB CORE
1658	Result Final Verified	Mary Ulmschneider	LAB CORE
1658	Reflex Action Occurred	Mary Ulmschneider	LAB CORE
1658	Result Entered	Mary Ulmschneider	LAB CORE
1658	Result Entered	Mary Ulmschneider	LAB CORE
1658	Test Transferred	Mary Ulmschneider	LAB CORE
1653	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1653	Allquot Created	Mary Ulmschneider	LAB SIC
1653	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1653	Result Entered	Mary Ulmschneider	LAB SIC
1652	Received	Mary Ulmschneider	LAB CORE
1652	Order Canceled On Instrument	Mary Ulmschneider	LAB CORE
1651	Collection Updated	Mary Ulmschneider	LAB CORE
1651	Order Sent To Instrument	Mary Ulmschneider	LAB CORE
1651	Specimen Created	Mary Ulmschneider	LAB CORE

171M-276T0002.2 LT BLUE-2.7mL - Lupus Inh Instrument ID:16561
 171M-276T0002.3 LT BLUE-2.7mL - Heparin Instrument ID:16562

- If there are additional containers on a specimen (.2, .3etc) you can click on those to expand the tracking for each:

171M-276T0002.2 LT BLUE-2.7mL - Lupus Inh Instrument ID:16561

Date	Detail	User	Location
10/03/2017			
1700	Result Filed To Chart	Service Account, Beaker Results Out	LAB SIC
1700	Charge Triggered	Mary Ulmschneider	LAB SIC
1700	Result Final Verified	Mary Ulmschneider	LAB SIC
1700	Test Canceled	Mary Ulmschneider	LAB SIC
1700	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1700	Result Entered	Mary Ulmschneider	LAB SIC
1700	Result Entered	Mary Ulmschneider	LAB SIC
1659	Charge Triggered	Mary Ulmschneider	LAB CORE
1659	Result Filed To Chart	Service Account, Beaker Results Out	LAB CORE
1658	Result Final Verified	Mary Ulmschneider	LAB CORE
1658	Reflex Action Occurred	Mary Ulmschneider	LAB CORE
1658	Result Entered	Mary Ulmschneider	LAB CORE
1658	Result Entered	Mary Ulmschneider	LAB CORE
1658	Test Transferred	Mary Ulmschneider	LAB CORE
1653	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1653	Reflex Action Occurred	Mary Ulmschneider	LAB SIC
1653	Result Entered	Mary Ulmschneider	LAB SIC
1652	Received	Mary Ulmschneider	LAB SIC
1652	Packing List Shipped	Mary Ulmschneider	LAB CORE
1652	Order Canceled On Instrument	Mary Ulmschneider	LAB CORE
1652	Packing List Locked	Mary Ulmschneider	LAB CORE
1652	Added To Packing List	Mary Ulmschneider	LAB CORE
1652	Received	Mary Ulmschneider	LAB CORE
1651	Collection Updated	Mary Ulmschneider	LAB CORE
1651	Order Sent To Instrument	Mary Ulmschneider	LAB CORE
1651	Specimen Created	Mary Ulmschneider	LAB CORE

171M-276T0002.3 LT BLUE-2.7mL - Heparin Instrument ID:16562

- Below the Tracking section you will see results for each test done on the specimen along with all the comments added, the ref ranges and the method (instrument) used to run the test. The green indicates which components were sent to the chart. On the test level it indicates the test has been verified. The red indicates there is an abnormal result and the arrow indicates whether

it was high or low. (If there were two exclamation points this would indicate a result was critical):

Lupus Inhibitor Evaluation

Res	Component	Value	Units	I	A	L	IE	R	Ref. Range	Method	Chart	PV
1	INTERPRETATION	Consistent With Lupus Like Inhibitor		!					Negative for Lupus Like Inhibitor...	EVOLUTION B	✓	
1	APTT PATIENT	63.7	sec	!-					26.0-36.0	EVOLUTION B	✓	
1	HEXAGONAL PL APTT	12.8	sec	!-					0-<8.0	EVOLUTION B	✓	
1	DWVT	37.8	sec							EVOLUTION B	✓	39.7
1	DWVC		sec							EVOLUTION B	✓	
1	HEX LA1	64.1	sec							EVOLUTION B	✓	
1	HEX LA2	51.3	sec							EVOLUTION B	✓	
1	STAT CHARGE									EVOLUTION B	✓	
1	DRVV SCREEN RATIO	0.94							0.00-1.20	EVOLUTION B	✓	10.00

Heparin, Either Std/Lmw

Res	Component	Value	Units	I	A	L	IE	R	Ref. Range	Method	Chart	PV
1	Heparin, STD LMW	0.22	U/mL							EVOLUTION A CORE	✓	
1	Heparin Type (Do Not Result)									EVOLUTION A CORE	✓	

- As you scroll through Specimen Inquiry you will see sections that display **Flags&Holds**, **Contacts**, **Order Details** as well as sections showing where, when & to whom results were sent. You can also create a report to fax on request from the **Create New Report** section.

Specimen Inquiry: 17IM-276T0002

CC Recipients for Orders on 17IM-276T0002 *Modify
LUPUS INHIBITOR EVALUATION, PLASMA
None Listed
HEPARIN, EITHER STANDARD / LMW, BLOOD
None Listed
APTT 1:1 MIX, PLASMA
None Listed

Results Sent Via Interface
Recipient
OHS LAB RES OUT Interface
OHS LAB RES OUT Interface
OHS LOINC RES OUT Interface

Create New Report
*Create requisition report
*Create specimen report

Printed Labels and Documents
Type
Specimen Label
OHSULab Zebra aliquot label
Specimen Label
OHSULab Zebra aliquot label

Charge Summary

Create New Report

Select recipients [1]

Print	Recipient	Type	Fax #	Device
<input type="checkbox"/>	Beaker, Leo	Patient		
<input checked="" type="checkbox"/>	Thomas Yackel, MD	Primary Care Physician		DHB7065P
<input type="checkbox"/>	PEACE HEALTH LABORATOR	Submitter		DHB7065P

Select All Clear All

Select orders for recipient Thomas Yackel, MD [2]

Print	Date/Time	Procedure	Status
<input checked="" type="checkbox"/>	10/03/2017 1651	LUPUS INHIBITOR EVALUATION REFLEX TO APTT MIX	Final result
<input checked="" type="checkbox"/>	10/03/2017 1651	HEPARIN, EITHER STANDARD / LMW, BLOOD	Final result

Select All Clear All Preview

Override printer: Print Cancel

- The **Printed Labels & Documents** section has details about where labels printed and when. The **Charge Summary** section shows a list of all the charges for each test on the specimen. You can click

on these charges to see more detail about them.

Printed Labels and Documents							
Type	Printer	Printed					
Specimen Label	DHBA318P LABELS	10/03/2017 1653 by Mary Ulmschneider				*Reprint	
OHSUlab Zebra aliquot label							
Specimen Label	DHBA318P LABELS	10/03/2017 1651 by Mary Ulmschneider				*Reprint	
OHSUlab Zebra aliquot label							

Charge Summary							
Lupus Inhibitor Evaluation 17IM-276T0002							
Proposed Charges	Qty	Billed Charges	Qty	Mod	Svc Date	Status	
85730 (CPT®) [HB-LAB APTT]	1	85730 (CPT®) [HB-LAB APTT]	1		10/03/2017	Filed to Resolute Hospital Billing	
85613 (CPT®) [HB-LAB DVVT]	1	85613 (CPT®) [HB-LAB DVVT]	1		10/03/2017	Filed to Resolute Hospital Billing	
85598 (CPT®) [HB-LAB HEXAGONAL PL APTT]	1	85598 (CPT®) [HB-LAB HEXAGONAL PL APTT]	1		10/03/2017	Filed to Resolute Hospital Billing	
Heparin, Either Std/Lmw 17IM-276T0002							
Proposed Charges	Qty	Billed Charges	Qty	Mod	Svc Date	Status	
85520 (CPT®) [HB-LAB HEPARIN, EITHER STD/LMW]	1	85520 (CPT®) [HB-LAB HEPARIN, EITHER STD/LMW]	1		10/03/2017	Filed to Resolute Hospital Billing	
APTT 1:1 Mix 17IM-276T0002							

9. At the end of the report are details about the collection of the specimen, showing what instructions appear at collection as well as any questions and the corresponding answers entered at collection.

Heparin, Either Std/Lmw	
Collection Question	Answer
Test name	
Collection Instructions	
Tube must be full	

APTT 1:1 Mix	
Collection Instructions	
Tube must be full	
Shipping Instructions	
Spin, FRZ 2, -20	

© 2017 Epic Systems Corporation. All rights reserved. PROPRIETARY INFORMATION - This item and its contents may not be accessed, used, modified, reproduced, performed, displayed, distributed or disclosed unless and only to the extent expressly authorized by an agreement with Epic. This item is a Commercial Item, as that term is defined at 48 C.F.R. Sec. 2.101. It contains trade secrets and commercial information that are confidential, privileged and exempt from disclosure under the Freedom of Information Act and prohibited from disclosure under the Trade Secrets Act. After Visit Summary, Analyst, App Orchard, ASAP, Beaker, BedTime, Bones, Break-the-Glass, Caboodle, Cadence, Canto, Care Everywhere, Charge Router, Chronicles, Clarity, Cogito ergo sum, Cohort, Colleague, Community Connect, Cupid, Epic, EpicCare, EpicCare Link, Epicenter, Epic Earth, EpicLink, EpicWeb, Good Better Best, Grand Central, Haiku, Happy Together, Healthy People, Healthy Planet, Hyperspace, Identity, IntraConnect, Kaleidoscope, Limerick, Lucy, MyChart, OpTime, OutReach, Patients Like Mine, Phoenix, Powered by Epic, Prelude, Radar, RedAlert, Resolute, Revenue Guardian, Rover, SmartForms, Sonnet, Stork, Tapestry, Trove, Welcome, Willow, Wisdom, and With the Patient at the Heart are registered trademarks, trademarks or service marks of Epic Systems Corporation in the United States of America and/or other countries. Other company, product and service names referenced herein may be trademarks or service marks of their respective owners. U.S. and international patents issued and pending.

This guide is based on Epic 2017.